

USER MANUAL

Tie Me Ghana

Hospital communication for Deaf and Hard of Hearing patients, in Ghanaian Sign Language.

Pilot release, September 2026

Built for the MTN Ghana Tekyerema Pa Hackathon 2026

Read this first

This app is not an interpreter and does not replace one. Where a qualified Ghanaian Sign Language interpreter is available and the conversation matters, use the interpreter.

This is a pilot. It has not been approved as a medical device and is not cleared for clinical use. Do not use it alone to take consent for a procedure, to deliver a serious diagnosis, or for any exchange where a misunderstanding carries real risk.

What this app does

Tie Me Ghana helps a clinician and a Deaf or Hard of Hearing patient understand one another during a hospital visit, without an interpreter in the room.

You type or speak. The patient sees your words as a video of a real person signing in Ghanaian Sign Language. They answer by tapping, typing or pointing, and their answer is read aloud to you so you can keep your hands and eyes on the patient.

How the signing works, and where it stops

Every sign is a recording of a real person, approved in advance by a fluent Ghanaian Sign Language consultant. The app never invents or animates a sign, so it can only show what has been filmed and approved.

The app refuses rather than guesses. When a word has no approved sign, it tells you instead of showing a different one. A near miss in a clinical sentence is more dangerous than a gap, because it looks like an answer.

Some sentences will be refused outright and will not be shown to the patient. When that happens, rephrase the sentence or find another way to ask.

What is in this manual

1. The opening screen, and starting a visit
2. Talking to a patient who reads
3. Talking to a patient who does not read
4. Emergency mode
5. Writing a prescription
6. Working with a poor connection
7. What the app stores
8. If something goes wrong


Before the first patient

- Open the app once on a good connection and let it finish setting up. It saves the sign videos onto the device so they play instantly afterwards, and shows its progress while it does.
- Tap **Install app** so it has its own icon and opens without a browser.
- Read the privacy policy and terms, linked at the bottom of the opening screen.

SECTION 1

The opening screen

This is what the app opens into, and what the next patient should always be handed. One question decides how the whole visit will work.

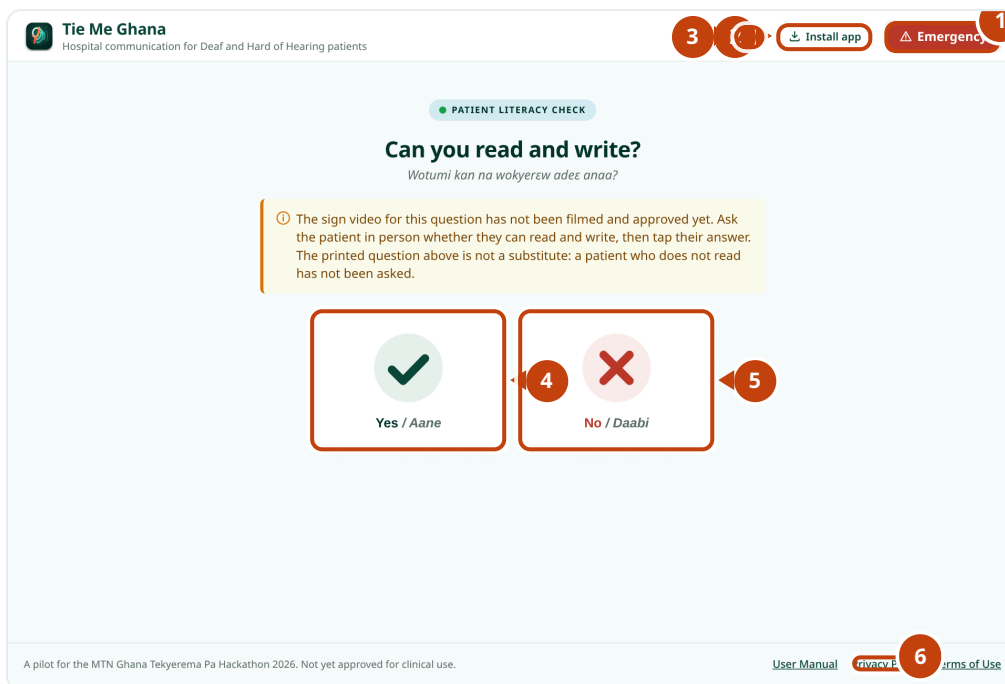

What each control does

1	Emergency	Opens Emergency mode straight away, without asking anything first. Use it for a patient who arrives after an accident or who cannot breathe. It works before you have asked whether they read.
2	Install app	Adds Tie Me Ghana to the device like any other app, with its own icon. Worth doing once on every device that will be used in the clinic: an installed copy opens faster and keeps its saved sign videos.
3	Connection	Green means the hospital system is reachable. If it turns grey the app keeps working with the sign videos already on the device, but new prescriptions cannot be issued until it returns.
4	Yes	The patient reads written language. Tapping this opens the full consultation screen, where you type freely and they can type back.
5	No	The patient does not read print. Tapping this opens Guided Interrogation, where every question is asked in sign video and answered yes or no.
6	Privacy and Terms	The full privacy policy and terms of use. Worth reading once before the app is used with real patients, and the page to send a hospital that asks what the app stores.

How to use it


1. Ask the patient, in person, whether they can read and write. The printed question on screen is not a substitute: a patient who does not read has not been asked.
2. Tap Yes or No on their behalf, or let them tap it.

3. The app opens the matching screen and stays there for the rest of the visit.

SECTION 2

Talking to a patient who reads

For a patient who reads. You type, they see both the words and the sign video, and they can type back.

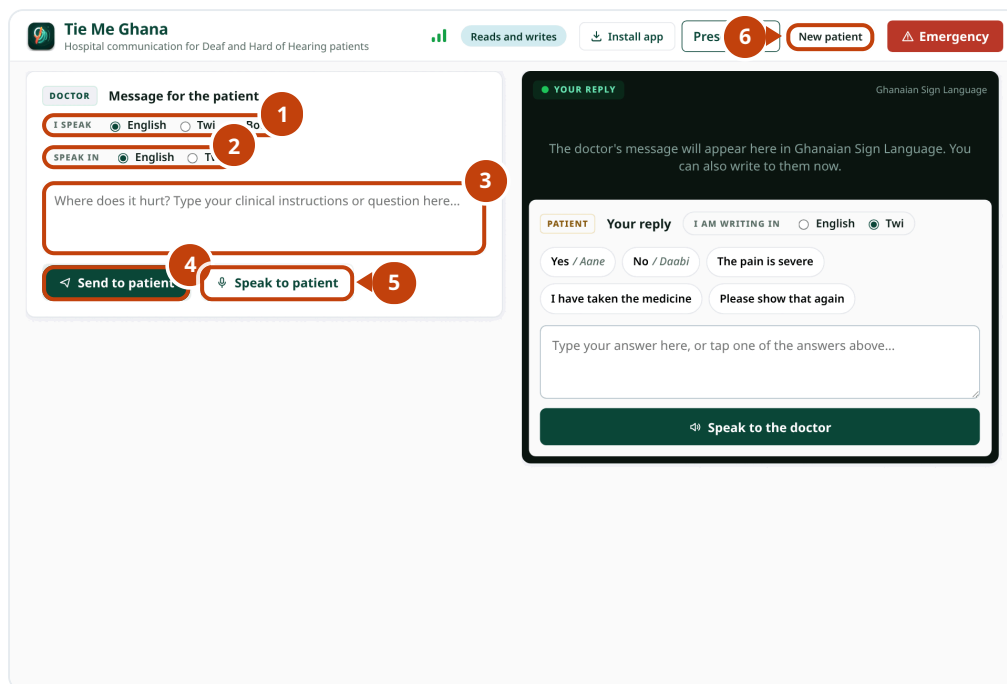

What each control does

1	I speak	The language you are about to use. Choose Both when your sentence mixes English and Twi, which is normal for medicine names that have no Twi word. Speaking aloud needs one language chosen, because speech recognition handles one at a time.
2	Speak in	The language the patient's answers are read aloud to you in. Set it once at the start of the visit to whichever language you understand.
3	The message box	Type what you want to say. Keep it to one clear instruction or question at a time. The app will show you the sign video before the patient sees it.
4	Send to patient	Renders your message as Ghanaian Sign Language and shows it to the patient. If a word has no approved sign, the app tells you instead of guessing.
5	Speak to patient	Dictate instead of typing. Tap once to start, tap again to stop and send. Typing always works and needs no permission, so use it whenever the room is noisy.
6	New patient	Ends this visit and erases the record of it from the device. Tap it before handing the device to the next patient. This is the control that stops one patient seeing another patient's consultation.


How to use it

1. Set I speak and Speak in once, at the start of the visit.
2. Type one clear instruction or question.
3. Tap Send to patient, then turn the device towards them.
4. Read their reply, or listen to it read aloud.

SECTION 3

Talking to a patient who does not read

For a patient who does not read print. Every question is asked in sign video and answered with a tap, so nothing depends on written language.

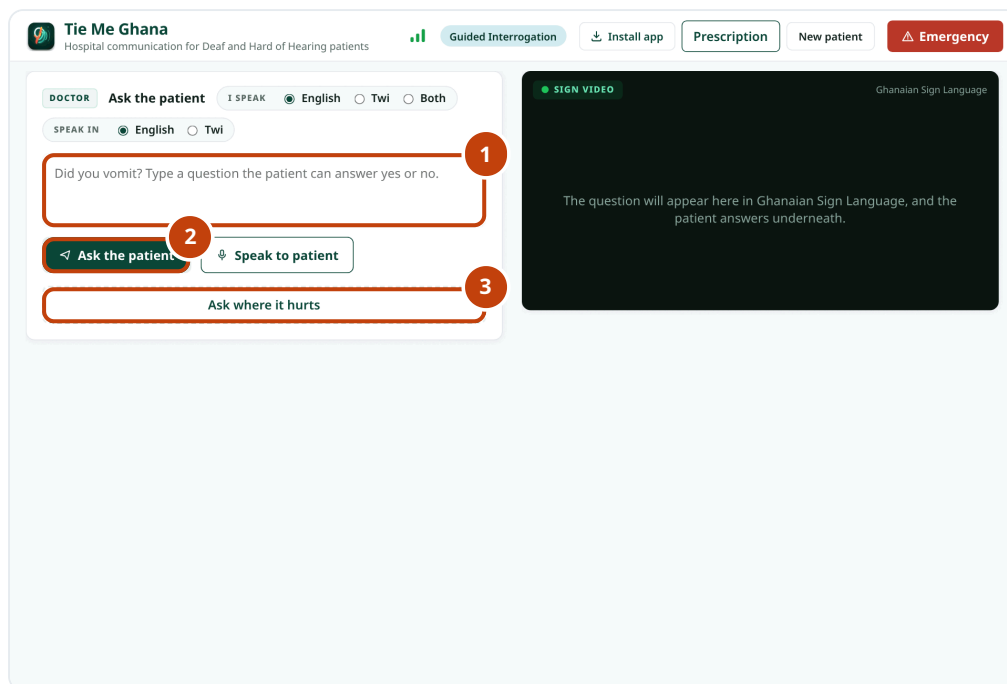

What each control does


1	The question	Ask one thing at a time. Guided Interrogation is for a patient who does not read, so every question becomes a sign video and is answered yes or no.
2	Send to patient	Shows the question as sign video. When the video ends, the Yes and No answers appear in the same place, so the patient does not have to look elsewhere.
3	Ask where it hurts	Asks the question and then offers the patient a set of body locations to choose from, instead of a yes or no answer.

How to use it

1. Type a question that can be answered yes or no.
2. Tap Send to patient. The question plays as sign video.
3. When the video ends, Yes and No appear where the video was. The patient taps one.
4. Use Ask where it hurts when you need a body location instead of yes or no.

Emergency mode

For an accident, a collapse, or any moment where there is no time to set anything up. Nothing here needs typing, reading, or a sign video that has been filmed.

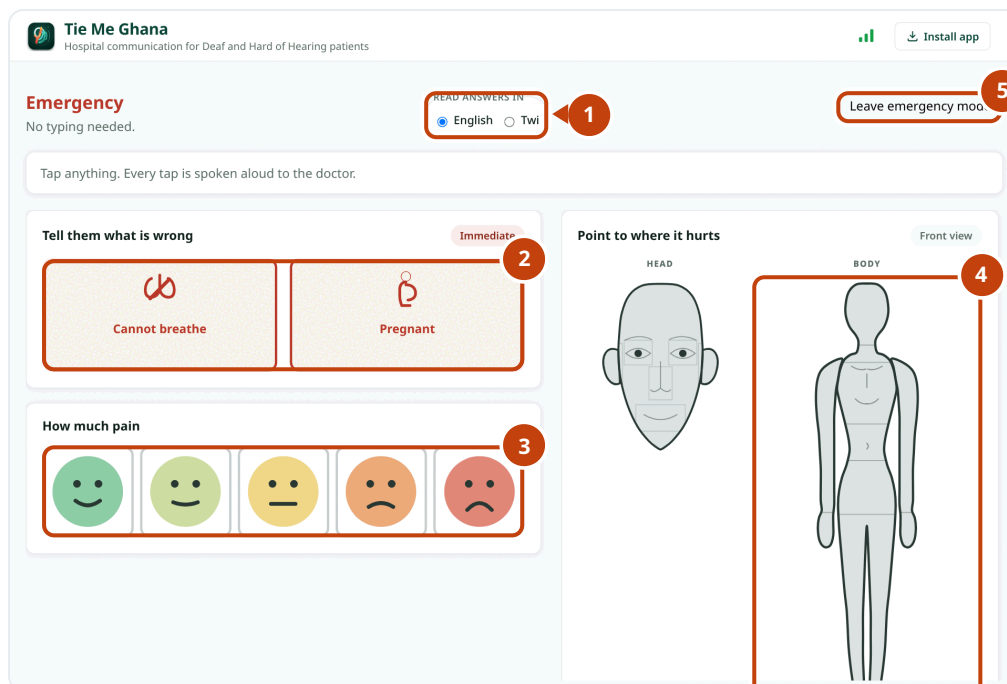

What each control does

1	Read answers in	The language every tap is read aloud in. Set it first, in the few seconds before you start, so what the patient says reaches you in a language you understand.
2	Critical alerts	One tap says the most urgent things a patient may need to tell you. Each is spoken aloud immediately. Where the sign video has not been filmed yet, the icon and the label are still shown and still work.
3	How much pain	Five faces, from comfortable to severe. The patient taps one and it is spoken aloud. This needs no sign video and no reading at all.
4	Point to where it hurts	The patient taps a part of the head or body and it is spoken aloud. Only the part they touched is highlighted, so there is no doubt about what they meant.
5	Leave emergency mode	Returns to the normal consultation. Emergency mode sits outside the visit, so leaving it does not lose anything you had already asked.

How to use it

1. Tap Emergency from any screen. It works before the literacy question has been asked.
2. Set Read answers in to the language you understand.
3. Hand the device to the patient. Every tap is spoken aloud to you immediately.
4. Tap Leave emergency mode when the immediate danger has passed.

Writing a prescription

The instructions the patient takes home, as a QR code they can replay in sign language on their own phone for as long as the course lasts.

Tie Me Ghana
Hospital communication for Deaf and Hard of Hearing patients

Reads and writes Install app New patient Emergency

Prescription

One row per medicine. The patient takes these home as sign language video.

Photograph the medicine
The patient sees this picture, then the dose in sign language. With a photo the name is optional.

Photograph the medicine itself, not a pharmacy label. A label may carry the patient's name, and this picture goes home with the QR code.

Medicine
Paracetamol

How much: 1 Of what: Tablet

WHEN
☐ Morning ☐ Afternoon ☐ Evening ☐ Night

How often: Twice a day Food: Any time For how long: days

THE PATIENT IS TOLD
one tablet, twice a day

Add another medicine Create QR code

What each control does

1	Photograph of the medicine	Take a picture of the packet. The patient recognises the box far more reliably than the name, and it means you do not have to type a drug name at all. Photograph the medicine, never the patient and never a label carrying a name.
2	How much	Half, one, two, and so on, with the unit beside it. These are chosen rather than typed so that every prescription the app can produce is one it can also sign.
3	When	Morning, afternoon, evening, night, and whether it is taken before, with or after food. Choosing specific times is clearer for the patient than a count such as twice a day.
4	Issue the prescription	Produces the QR code the patient takes home. Scanning it replays the whole list in sign language, on their own phone, with no account and no app to install.

How to use it

1. Photograph each medicine, so the patient recognises the packet rather than a name.
2. Choose how much, and when.
3. Tap Issue the prescription.
4. Show the patient the QR code so they can scan it, and print it if they want a paper copy.

Working with a poor connection

Hospital connections drop. The app is built for that, and most of it keeps working.

WHAT YOU ARE DOING	WITH NO CONNECTION
Asking a question the device has already shown	Works. The sign video is saved on the device.
Asking something new	Needs a connection to look up the signs.
Emergency mode: pain scale and body map	Works. Both are drawings and need nothing from the server.
Reading an answer aloud	Needs a connection.
Issuing a prescription	Needs a connection.
A patient replaying a prescription at home	Works once they have opened it once on a connection.

The indicator in the top bar

Green means the hospital system is reachable. Grey means it is not. The app does not stop when it turns grey: it keeps using what is already on the device and tells you which actions are unavailable rather than failing silently.

On the day of a demonstration or a clinic session, open the app a minute beforehand. The pilot runs on free hosting that puts the server to sleep after about fifteen minutes of no use, and waking it takes up to a minute.

What the app stores

The short version. The full privacy policy is linked at the bottom of the opening screen.

INFORMATION	WHERE IT GOES
The patient's name, age, number or diagnosis	Never collected. There is nowhere in the app to enter them.
The record of the visit	The device only. It is never sent anywhere, and it is erased when you tap New patient.
Prescriptions	Stored on the server as medicine names, photographs and dosages, with nothing saying whose they are.
Photographs of medicines	Stored with the hidden camera information removed, including the location the photograph was taken.

Two habits that matter.

Tap **New patient** before handing the device on. That is what erases the previous consultation, and these consultations cover pregnancy, sexually transmitted infections and HIV status.

Photograph the medicine, never the patient, and never a label carrying somebody's name.

The prescription QR code

Anybody who scans the code or holds the link can see that list of medicines. It carries no name, so it says what somebody is taking without saying who, but treat a printed code the way you would treat any prescription slip.

SECTION 8

If something goes wrong

The problems most likely to come up, and what to do about each.

WHAT YOU SEE	WHAT IT MEANS, AND WHAT TO DO
The app says a word has no sign	That sign has not been filmed and approved yet. Rephrase using simpler words, or write the word down if the patient reads.
A question is refused and will not send	The sentence contains something that cannot be signed safely. The app will not show a half correct clinical sentence. Rephrase it.
The video area says it is getting ready	The sign video is still downloading. It plays as soon as it arrives. If this happens often, open the app on a better connection once and let it finish setting up.
The connection indicator is grey	The hospital system is unreachable. Saved sign videos, the pain scale and the body map still work. New questions and prescriptions do not.
The microphone will not start	Either the browser refused permission, or Both is selected as your language. Speech recognition handles one language at a time, so choose English or Twi, or type the message.
The first request after a quiet period is slow	The pilot server sleeps when unused and takes up to a minute to wake. Open the app a minute before a clinic session.
Captions are not translated and spoken answers are silent	Expected in this pilot. Translation is switched off and the app says so on screen rather than pretending to work.

Getting help

Write to bestdrtrick@gmail.com. Include what you were doing, what you expected, and what happened instead. If it concerns a prescription, include its reference, which is the code under the QR image.